Troubleshooting **Dépannage**

When a malfunction or incorrect temperature measurement occurs, an error message will appear. The thermometer will beep four times, then automatically turn off.

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Technical Specifications **Spécifications techniques**

Measurement Range: 35°C to 42°C (95°F to 107.6°F).
Measurement Accuracy: Technical: +/- 0.1°C (0.2°F)
Operation Environment: 10°C to 40°C (50°F to 104°F),
 Humidity: 15-95% (non condensing) (complies with ASTM E 1112-86)
Storage Environment: -20° C to 50° C (-4°F to 122°F),
 Humidity: 15-95% (non condensing)

Test Protocol:
1 Install the oral or rectal probe to the handle holder.
2 Put the oral or rectal probe into water bath.
3 Wait for 3 seconds, and then turn on the thermometer by pressing the O/I button (the thermometer has been in standard mode for measurement). When thermometer determined the peak value of temperature, the backlight turns on with two beeps. Temperature measurement completed and then records the value of temp.
4 For the underarm probe, it shall be inserted to water bath for 60 seconds prior and then turns on the thermometer for measurement. When the thermometer determined the peak value of temp., the backlight turns on with two beeps. Measurement completed and records the value of temperature.

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Dorel Juvenile Group, Inc. warrants to the original purchaser that this product (Complete Family Thermometer), is free from material and workmanship defects when used under normal conditions for a period of one (1) year from the date of purchase. Should the product contain defects in material or workmanship Dorel Juvenile Group, Inc. will repair or replace the product, at our option, free of charge. Purchaser will be responsible for all costs associated with packaging and shipping the product to Dorel Juvenile Group Consumer Relations Department at the address noted on the front of this document and all other freight or insurance costs associated with the return. Dorel Juvenile Group will bear the cost of shipping the repaired or replaced product to the purchaser. Product should be returned in its original package accompanied by a proof of purchase, either a sales receipt or other proof that the product is within the warranty period. This warranty is void if the owner repairs or modifies the product or the product has been damaged as a result of misuse.

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This warranty excludes any liability other than that expressly stated previously, including but not limited to, any incidental or inconsequential damages.

©2008 Dorel Juvenile Group. All Rights Reserved.
 www.djgusa.com (800) 544-1108 www.safety1st.com
 Styles and colors may vary. Les styles et les couleurs peuvent varier.
 Distributed by (distribué par) Dorel Juvenile Group, Inc., 2525 State St., Columbus, IN 47201-7494
 Dorel Distribution Canada, 873 Hodge, St. Laurent, QC H4N 2B1
 11/20/08 4358-4764

49546 User Guide
Complete Family Thermometer

! WARNING:
BEFORE USING DIGITAL THERMOMETER, PLEASE READ CAREFULLY.
 • This is **NOT** an ear thermometer. Take temperature measurements Underarm, Orally or Rectally.
 • **NEVER** allow a child to play with these products. They are not toys. They are sensitive electronic devices and may be damaged if not used properly under parental supervision.
 • **ALWAYS** supervise your child when using thermometer. **NEVER** allow your child to use it alone or without close parental supervision.
 • Keep battery out of child's reach. Swallowing battery could be fatal. Battery should not be charged or placed into extreme heat as it may explode.
 • Used batteries should be disposed of properly. Do not dispose of batteries in fire, recharge, put in backwards, disassemble, mix with used or other battery types – they may explode or leak and cause injury.
 • Avoid sudden shock and dust.

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
Important Information **Information importante**

• Read these instructions carefully to ensure accurate temperature readings. Temperature readings are affected by many factors including exercise, drinking hot or cold beverages prior to taking a reading as well as technique.
 • High, prolonged fever requires medical attention especially for young children. Please contact your physician.
 • Do not allow children to take their temperatures unattended. This thermometer is not a toy and should be kept out of a child's reach when not in use.
 • Do not allow children to walk or run during temperature taking. The probe could present a hazard.
 • Oral, underarm, and rectal temperatures vary by approximately 0.5°C (1.0°F).
 • Do not return this product to the place of purchase. If any parts are missing, email consumer@djgusa.com, call Consumer Relations at (800) 544-1108, or fax at (800) 207-8182. You can also visit our website at www.safety1st.com. Have the model number (49546) and date code (manufacture date located on back of thermometer) ready.

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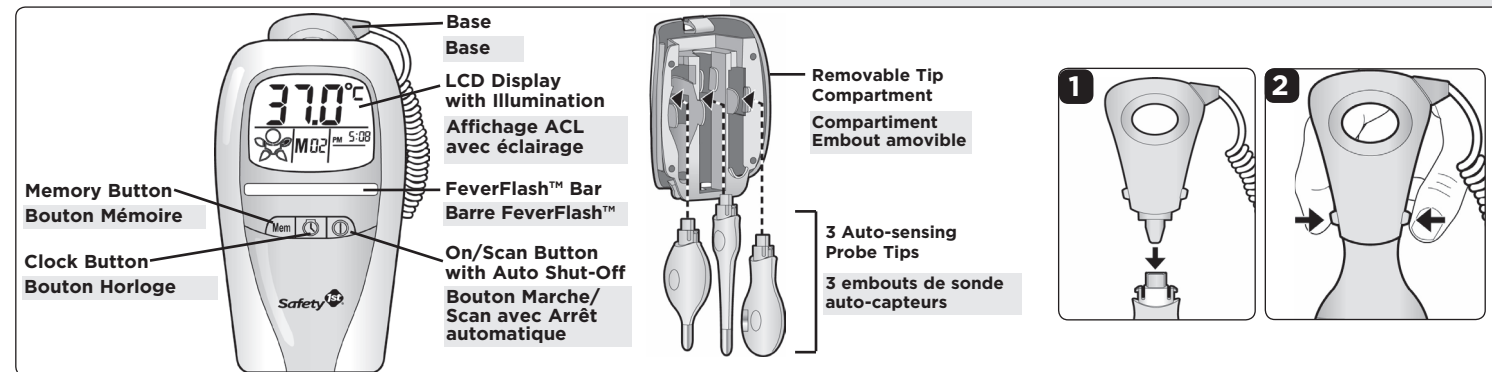1 To attach a tip, select oral, rectal, or underarm tip. Align probe tip base to plug tip (Figure 1). Tip fits in one direction only. The mode is automatically sensed and will display accordingly.
2 To remove tip, press buttons on both sides and pull tip off (Figure 2).

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Temperature Measurements Range

Many people believe that normal body temperature is 37.0°C (98.6°F). However, medical research has shown that temperature varies from person to person. Factors affecting body temperature include age, gender, physical activity, environmental temperature, clothing, medications, and time of day. Oral, rectal, and underarm temperatures vary by approximately 0.5°C/ 1.0°F.

| Method | Average Normal Temperature Range | Fever* |
|----------|--------------------------------------|-----------------------|
| Oral | 36.1°C to 37.8°C (97.0°F to 100.0°F) | Over 37.8°C (100.0°F) |
| Rectal | 36.3°C to 38°C (97.3°F to 100.4°F) | Over 38°C (100.4°F) |
| Underarm | 35.2°C to 37.2°C (95.4°F to 99.0°F) | Over 37.2°C (99.0°F) |

*Fever definitions according to the American Medical Association

Set Clock


NOTE: Set the clock at time of first use, after replacing battery, or any time the clock button is accidentally pressed. **To set clock:** Press clock button . Hours will blink. Press memory button **Mem** until correct hour displays, then press clock button again. Minutes will blink. Press memory button until correct minutes display, then press clock button again. "AM" will blink. Press memory button to scroll between AM or PM. When time is correct, press the on/scan button .

Changing from Fahrenheit to Celsius

To switch the display between °F and °C and vice-versa: While the thermometer is off, press and hold the On/Scan button for about 5 seconds until you see the °F or °C. Release button and quickly press it again. Unit has been reset.

FeverFlash™

The FeverFlash™ occurs during temperature measurement when the temperature is high. (See fever chart included on package or under "Temperature Measurement Range"). If your child has an elevated temperature, the red FeverFlash™ bar will light up and multiple beeps will sound. **High or prolonged fever requires medical attention, especially for young children. Contact your physician.**

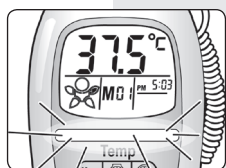Memory Function

The thermometer memory function stores the last 10 measurements along with the time. The last measurement will appear for approximately 3 seconds after turning the thermometer on and after it runs its built-in test. This feature allows you to compare your most recent measurement to the prior ones. **To scroll through the last 10 measurements:** With the thermometer off, press and release the memory button **Mem**. The screen will display the number of the reading (such as **MB 1**) along with that measurement and the time that measurement was taken. Press the memory button to scroll to reading number 2. Press the memory button again, and it will scroll up to number 3, and so on. When 10 spaces are filled the oldest measurement will be the next to be deleted.

Operating Thermometer

Getting the Thermometer Ready:

1 Two Ways to Turn On: Press and release the On/Scan button , or remove the probe tip base from the top of the unit. One beep occurs and all segments are displayed with back light for two seconds (**Figure 1**). The previous measurement is displayed for three seconds (**Figure 2**). Backlight shuts off and "- -" flashes until step 2 is complete.

2 Select Mode: Open the door on the back of the base (see "removable tip compartment" under section "Know Your Unit"). Select oral, rectal, or underarm tip and plug into probe tip base (**Figure 3**). The mode is automatically sensed and will display accordingly (**Figure 3**). **NOTE:** Tip fits in one direction only.

3 When you hear two beeps, the thermometer is ready to measure temperature.

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Operating Thermometer (continued)

Choosing a Method:

ORAL METHOD—For children 4 years and older.

1 Place the probe tip under the tongue in locations indicated by (**Figure 1**). Instruct children not to bite down on the probe.

UNDERARM METHOD

1 Wipe underarm with a dry towel. Probe must touch skin.

2 Place probe tip at center of underarm (**Figure 2**) and bring arm down against the body.

RECTAL METHOD—Commonly used for babies, young children or when it is difficult to take an oral or underarm temperature.

1 Lubricate probe tip with a water-soluble lubricant or petroleum jelly.

2 For easier insertion, position baby so his back and legs form a 90-degree (L-shaped) angle. Try placing child on stomach on your lap with his legs down. Gently insert the probe tip no more than 1/2 inch (1.3 cm) into the rectum. Do not force the tip into the rectum if resistance is encountered.

Taking the Temperature:

1 Once you have positioned the thermometer, hold it steady while the temperature is obtained. One long beep will be heard to indicate that measurement is complete. Back light will illuminate and "clean me!" icon will blink. Temperature will be displayed for about one minute without back light. **FEVER ALERT NOTE:** If patient has an elevated temperature, reading red FeverFlash™ bar light up for 8 seconds and multiple beeps in ten seconds will sound.

2 Press the On/Scan button to turn the thermometer off. If left on, automatic shut off will occur in approximately one minute.

NOTE: Wait at least one minute before taking another measurement. Clean the thermometer (see "Care and Maintenance" section).

Care and Maintenance

• Clean the probe and tip with an alcohol swab before and after each measurement.

• Use a soft, dry cloth to clean thermometer body. Never use abrasive cleaning agents, thinners or benzene for cleaning. Do not scratch the surface of the display. Do not expose the thermometer to extreme temperatures, humidity, direct sunlight or shock. Store at room temperature.

To Replace Battery

When the "Low Battery" symbol appears (**Figure 1**), the battery is exhausted and needs replacing.

1 Use a small screwdriver to loosen screw. Slide cover off (**Figure 2**).

2 Insert 2 AAA alkaline batteries (**Figure 2**). Position batteries correctly with the positive "+" side as indicated in the battery compartment.

3 Slide the battery cover on and secure the screw.

NOTE: Properly dispose of the batteries, keeping them away from small children and heat (**see Warning section**).

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