Hospital's Choice

## Complete Healthcare Kit 62403 & 49373

### Digital Thermometer with Case

Model 10430

## WARNING:

**BEFORE USING DIGITAL THERMOMETER, PLEASE READ CAREFULLY.**

- **POTENTIAL CHOKING HAZARD.** Do not leave child unattended when a thermometer cover is in use.
- Thermometer covers are not recommended for oral use.
- Keep thermometer covers out of reach of children.
- Thermometer covers are intended for single use only.
- This is **NOT** an ear thermometer. Take temperature measurements Underarm, Orally or Rectally.
- **NEVER** allow a child to play with these products. They are not toys. They are sensitive electronic devices and may be damaged if not used properly under parental supervision.
- **ALWAYS** supervise your child when using thermometer and thermometer covers. **NEVER** allow your child to use it alone or without close parental supervision.
- Keep battery out of child's reach. Swallowing battery could be fatal. Battery should not be charged or placed into extreme heat as it may explode.
- Used batteries should be disposed of properly. Do not dispose of batteries in fire, recharge, put in backwards, disassemble, mix with used or other battery types – they may explode or leak and cause injury.

The Safety 1<sup>st</sup> Digital Thermometer takes a clear and accurate temperature without the worries of a glass/mercury device. Effective for oral, rectal or underarm use.

### BEFORE USING THE DIGITAL THERMOMETER, PLEASE READ CAREFULLY:

- This thermometer is not a toy and should be kept out of a child's reach when not in use.
- Keep thermometer away from hot or boiling water, chemicals, and direct sunlight.
- Do not chew on or bend the thermometer.
- Make sure the thermometer has been at room temperature for 15 minutes before use.

### Important Information:

- High, prolonged fever requires medical attention especially for young children. Please contact your physician.
- Do not allow children to take their temperatures unattended.
- Read the enclosed instructions carefully to ensure accurate temperature readings. Temperature readings and measurement times are affected by many factors including exercise, drinking hot or cold beverages prior to taking a reading as well as technique.
- Do not allow children to walk or run during temperature taking. The probe could present a hazard.
- Oral, underarm, and rectal temperatures vary by approximately 1°F (0.5°C).

### Special Product Features:

- Fast 60-second rectal and oral readings! 90-Second underarm readings.
- Stops flashing when reaches highest temperature.
- Automatic switch off after approximately 10 minutes- for longer battery life.
- Digital display, for temperature changes as minor as 0.2°F (0.1°C).
- Water resistant for easy cleaning.
- Long life, non-replaceable, mercury free battery (approximately 6000 temperatures).
- Low battery indicator “▼”.

### How To Care For This Product:

- Use of thermometer covers is recommended to help avoid spreading germs and to keep the thermometer clean. The thermometer should be cleaned after each use even when thermometer covers are used.
- Clean the thermometer tip and probe using isopropyl alcohol or soap and water.
- Store the unit in the protective case while not in use.
- Do not store or keep in direct sunlight, or in high temperature areas.
- Do not drop the thermometer, as this may ruin the unit.

### What Is A “Normal” Temperature?:

Temperature readings vary in different individuals. A normal temperature refers to each individual's average temperature while they are well. Normal average temperatures range from 97°F (36.1°C) to 100°F (37.7°C) when taken orally. From an oral temperature, you may calculate a rectal equivalent temperature by adding 1°F (0.5°C) while an underarm temperature requires subtracting 1°F (0.5°C). For example, if your oral temperature is 98.6°F (37°C) this is the same as a rectal temperature of 99.6°F (37.5°C) or an underarm temperature of 97.6°F (36.4°C).

### How To Take Temperature:

- 1 Select rectal, oral or underarm as your measurement site based on the individual requiring a temperature measurement.
- 2 Push the on-off button located on the cap. The thermometer display will read: “188.8 °F”. This is a function check indicating the unit is working properly. The thermometer will read L°F (or L°C) and begin measuring temperature when exposed to body heat.
- 3 The °F (or °C) symbol will continuously flash until completion of the temperature taking process. Please note that the “▼” indicates a low battery. An “ERR” message in the window indicates your unit is not functioning properly. Do not expect an accurate temperature reading if this low battery indicator is present.
- 4 When peak temperature is reached, the °F (or °C) symbol will stop flashing.
- 5 Read and record the temperature and time for reference.
- 6 Push the on/off button to turn the unit off. If you do not turn the unit off, it will shut off automatically in approximately 10 minutes.


### Rectal Use:

- Use of thermometer cover is recommended for rectal use.
- Clean the thermometer as recommended.
- Commonly used for babies, young children or when it is difficult to take an oral or underarm temperature.
- Lubricate the thermometer cover with a water-soluble lubricant or petroleum jelly.
- Gently insert the probe tip no more than 1/2 inch or 1.3 cm into the rectum. Do not force the tip into the rectum if resistance is encountered.
- The peak temperature should be reached in approximately 60 seconds.

**NOTE:** To purchase thermometer covers call Consumer Relations at (800) 544-1108, and ask for part number 49520.


### Oral Use: (For children 4 years and older)

- Do not drink hot or cold fluids, exercise or perform other activities that will raise or lower temperature readings when compared to your normal, average temperature. The mouth should remain closed up to 5 minutes before attempting a reading.
- Place the probe tip well under the tongue. (See diagram) as indicated by the “☺”. Instruct children not to bite down on the probe.
- The peak temperature should be reached in approximately 60 seconds. Opening the mouth could result in a longer time for a reading.
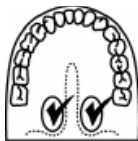
### Underarm (Axillary) Use:

- Wipe underarm with a dry towel.
- Place probe tip in armpit and bring arm down against the body. The thermometer should be pointing in the same direction as your body if possible.
- The time required for an accurate underarm temperature may be 90 seconds or longer depending on placement of the probe and closure of the arm against the body.
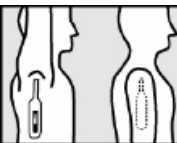
### Battery Information:

When the “▼” appears in the lower right of the display, the battery is exhausted. Dispose of the thermometer as the battery is non-replaceable. Be aware of any local regulations concerning disposal of this product. The battery is well sealed inside the thermometer to reduce the risk of being swallowed by children.

**Measurement Range:** 89.6°F to 111.2°F / 32°C to 44°C.

**Measurement Accuracy:** 0.2°F / 0.1°C between 93.2°F to 107.6°F / 34°C to 42°C at an ambient temperature of 64.4°F to 82.4°F / 18°C to 28°C.

°F	°C
96	35.6
97	36.1
98	36.7
99	37.2
100	37.8
101	38.3
102	38.9
103	39.4
104	40.0
105	40.6

## One Year Limited Warranty

Dorel Juvenile Group, Inc. warrants to the original purchaser that this product (Digital Thermometer), is free from material and workmanship defects when used under normal conditions for a period of one (1) year from the date of purchase. Should the product contain defects in material or workmanship Dorel Juvenile Group, Inc., will repair or replace the product, at our option, free of charge. Purchaser will be responsible for all costs associated with packaging and shipping the product to Dorel Juvenile Group Consumer Relations Department at the address noted on front page and all other freight or insurance costs associated with the return. Dorel Juvenile Group will bear the cost of shipping the repaired or replaced product to the purchaser. Product should be returned in its original package accompanied by a proof of purchase, either a sales receipt or other proof that the product is within the warranty period. This warranty is void if the owner repairs or modifies the product or the product has been damaged as a result of misuse.

This warranty excludes any liability other than that expressly stated previously, including but not limited to, any incidental or inconsequential damages.

## Healthcare Essentials

### Using Medicine Spoon & Dropper

- **Child should sit upright when taking medicine. Hold baby upright when administering medicine.**
- **Administer medicine slowly enough so your child may swallow comfortably. Pause and repeat if necessary.**
- **Wash thoroughly before first use and after each use.**

**CAUTION:** KEEP ALL MEDICATIONS OUT OF CHILD'S REACH. WE RECOMMEND THAT YOU USE A CHILD RESISTANT TOP TO SEAL THE MEDICINE BOTTLE AND STORE IT IN A LOCKED CABINET OR OTHER LOCATION THAT IS NOT ACCESSIBLE TO A CHILD.

### Medicine Spoon

Giving medicine to your child is easy with the Safety 1<sup>st</sup> Medicine Spoon. It offers an easy grip handle, smooth rounded edges and easy to read dosage marks.

### Use and Care:

- Hold Medicine Spoon straight up and pour medicine to prescribed dosage.
- Make sure that the child is sitting in an upright position before giving medicine.

### To Clean:

Wash with warm, soapy water before first use and after each use. Rinse thoroughly. This spoon is also top rack dishwasher safe.

### Medicine Dropper

Safety 1<sup>st</sup> understands that giving a child medicine takes special attention. Our Medicine Dropper is specially designed to provide accurate dosage readings. The convenient spill guard plug helps reduce messy spills at home or away.

### Use and Care:

- Remove spill guard from dropper tip and squeeze air out of dropper.
- Insert dropper into the medicine bottle and slowly release pressure on the bulb until the medicine reaches the desired dosage marking.
- Remove dropper from the medicine bottle and tip upside down while securing the spill guard over the dropper tip.
- Administer medicine while child is in an upright sitting position. Remove spill guard and point tip of dropper to the inside of baby's cheek and gently squeeze the bulb.

### To Clean:

Wash with warm, soapy water before first use and after each use. Rinse thoroughly.

### Nasal Aspirator

Congestion due to colds is uncomfortable for baby. The Safety 1<sup>st</sup> Nasal Aspirator will gently relieve baby's discomfort by removing excess mucus or discharge to help baby breath easier.

### Use and Care:

- Hold baby upright.
- Squeeze bulb of Aspirator to empty out the air.
- Gently place tip of Aspirator into baby's nose so that it points up or down, but not directly towards the center of the nostril.
- Slowly release the bulb and the suction will draw excess mucus from the baby's nasal passage.
- Repeat in other nostril.

### To Clean:

Clean before first use and after each use. Gently squeeze air from Aspirator. Rinse the bulb in warm, soapy water. Squeeze and release bulb several times to thoroughly wash interior of Aspirator. Rinse in the same fashion with warm, clean water. Make certain that all of the water is drained during drying.

**CAUTION:** NEVER push the aspirator too far into the nostril and always dry thoroughly after each use. NEVER use in baby's mouth or ears. When storing, keep away from heat and bright light.

### Nail Clippers

Baby's nails grow very fast and can become quite sharp. Because of this, it is very important to keep nails short and smooth so babies won't scratch themselves. The Safety 1<sup>st</sup> Nail Clippers are specially designed with curved edges for baby's little nails.

### Use and Care:

- Sterilize by wiping edges with rubbing alcohol.
  - Clip nails straight across. Do not try to cut rounded edges.
- NOTE:** It is best to trim after a bath when nails are soft, or at night when baby is sleeping.

**CAUTION:** Be careful not to trim too close to baby's sensitive skin.

## CAUTION:

THESE PRODUCTS ARE NOT TOYS. Do not allow children to play with them. To prevent accidental injury, your child should never be left unattended. Products should only be used with adult supervision. When not in use, keep out of reach of children.

©2008 Dorel Juvenile Group. All Rights Reserved.  
www.djgusa.com (800) 544-1108  
www.safety1st.com  
Made in CHINA. Fabriqué en CHINE.  
Styles and colors may vary. Les styles et les couleurs peuvent varier.  
Distributed by (distribué par ) Dorel Juvenile Group, Inc.,  
2525 State St., Columbus, IN 47201-7494  
Dorel Distribution Canada, 873 Hodge, St. Laurent, QC H4N 2B100  
04/04/08 4358-4562

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